



## CareGard Quick Reference

Press the Volume button repeatedly until the volume level is comfortable.

Power/Home/Away lights up when AC Power is on.

Press it to turn Activity Check and Medication Reminders off when you are Away and back on when you are Home. The button flashes while you are Away.

Press the CALL button or your personal help button to signal for help if you having an emergency.

If the Reset/Status button is flashing, press it and the panel will tell you the current status. If there is a problem, consult your User Guide.

For more detailed information on the operation of the CareGard system refer to the CareGard User Guide.

## Underside of CareGard

OFF  
ON


Plug the large end of the included phone cord into the **LINE** jack on the panel.  
Plug the small end into a phone jack in your home.  
Plug any phones or answering machines into the **PHONE** jack.

**SETTING THE CLOCK**  
Step 1: Enter your access code (default 1234) using the red numbered keys.

**On/Off Switch.**

**SETTING THE CLOCK**  
Step 2: Use the Clock Set, Hours and Minutes buttons to set the time.  
Follow the voice instructions.

Use the **MENU** button to access phone and sensor tests and program medication reminders.

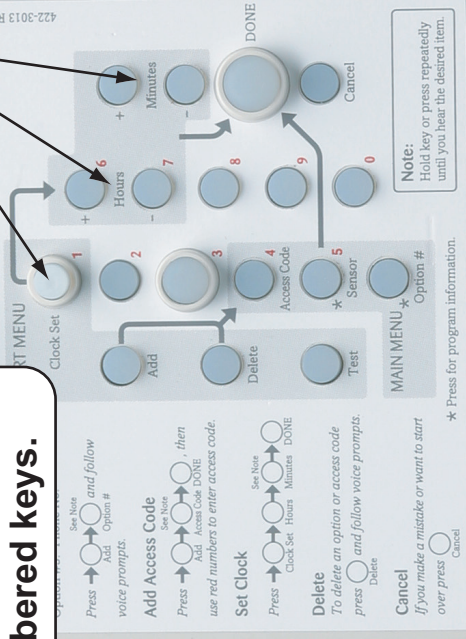

CareGard with cover removed



GE Interlogix

For more detailed information on the operation of the CareGard system please refer to the CareGard User Guide.